


TOE - Toe Separators

EN	TOE SEPARATOR
DE	ZEHENSEPARATOR
FR	SÉPARATEUR POUR ORTEILS
ES	SEPARADOR DE DEDOS DEL PIE
IT	SEPARATORE DI DITA DEI PIEDI
NL	TEENSCHIEDER
DK	TÅSEPARATOR
FI	VARPAANEROTIN
SE	TÅSEPARATOR
NO	SKILLER, TÅ
PL	SEPARATOR PALCÓW U STÓP
GR	ΔΙΑΧΩΡΙΣΤΗΣ ΔΑΚΤΥΛΩΝ ΠΟΔΙΟΥ

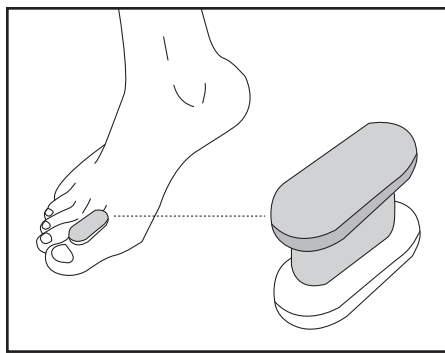

ENGLISH

Product Description

Silicone toe separator, reduces pain by preventing friction.

Usage Information (See diagram)

The Toe Separator prevents the first toe from overlapping the second, by being placed between the toe.

Indications for Use

- Prevents the first toe from overlapping the second which helps prevent deformity in the first metatarsophalangeal joint.
- Ideal for people with bunions or calluses.
- Reduces pain by preventing friction.

Contraindications

- Do not use over open wounds.
- Any known allergy and/or hypersensitivity to any material listed in product composition.
- Caution should be used if prescribed to patients with diabetes, vascular deficiency and neuropathy.

Warnings and Precautions

We recommend the initial fitting of this device be conducted by a suitably qualified healthcare professional who will advise the period of use.

- Carefully read all instructions and warnings prior to use.
- Follow all instructions to ensure proper performance of the device.
- Do not use if liniments, ointments, gels, creams or any other substances have been applied to the affected area.
- Do not re-use for another patient, doing so risks cross infection and can compromise product integrity.
- Perform regular skin and circulation checks, especially patients with diabetes, vascular deficiencies and neurological conditions.
- Should any adverse reactions occur, please cease use and contact your healthcare professional or provider of this product.
- The durability of the device may be compromised by certain factors, e.g., objects with sharp edges or damage to the hook and loop fasteners. To prevent this the hook and loop fasteners should always be fastened when the device is not being worn or when being washed.

Washing and Care Instructions

1. Hand wash in cool water using mild detergent.
2. Rinse and dry flat.
3. Do not spin or tumble dry.
4. Do not use fabric conditioner at any time.
5. Check for any degradation that may affect use.

Product Family Composition

Silicone.

Storage and Transport Conditions

Store in a cool, dry place out of direct sunlight, in the original packaging.

Recycling and Disposal

Packaging and constituent parts should be recycled or safely disposed in accordance with local or national laws.

Serious Incident

Report any serious incident to the manufacturer and the competent authority of the EU Member state, or country in which you reside.

MD MEDICAL DEVICE

1 SINGLE PATIENT - MULTIPLE USE

! CAUTION



Promedics Orthopaedics Ltd.
Port Glasgow Industrial Estate, Gareloch Road, Port Glasgow,
Scotland, UK, PA14 5XH
Tel: +44 (0) 1475 746 400 | Email: enquiries@promedics.co.uk
www.promedics.co.uk

CE **EC** **REP** **EMERGO EUROPE**
Westervoortsedijk 60,
6827 AT Arnhem
The Netherlands