


WTRS - Wrist Thumb Restriction Brace

EN	WRIST/HAND ORTHOSIS
DE	HANDGELENK-/HANDORTHESE
FR	ORTHÈSE DU POIGNET/DE LA MAIN
ES	ORTESIS DE MUÑECA/MANO
IT	ORTESI PER POLSO/MANO
NL	POLS-HANDORTHESE
DK	HÅNDELED-/HÅNDORTOSE
FI	RANNE-KÄMMENORTOOSI
SE	HÅNDELED- OCH HANDORTOS
NO	ORTOSE, HÅNDELEDD/HÅND
PL	ORTEZA NADGARSTKA/DŁONI
GR	ΝΑΡΘΗΚΑΣ ΚΑΡΠΟΥ/ΥΧΕΡΙΟΥ

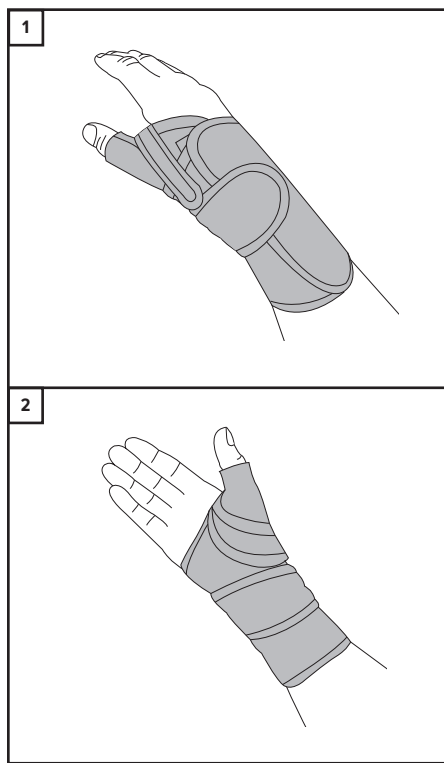

ENGLISH

Product Description

A ventilated wrist thumb brace with additional wrist compression wrap and thumb CMC abduction strap.

Usage Information (See diagram)

- Open brace and slide thumb into thumb sleeve and fasten hook/loop closures.
- Locate wrist strap and wrap around the wrist to give desired level of compression- there is an opening on the strap to aid grip for patients with limited dexterity.
- Pass the strap through the thumb web space and fasten to the location point on the back of the hand to resist medial subluxation.

Indications for Use

- Combined wrist and CMC joint arthritis.
- CMC medial subluxation.
- Tendonitis.

Contraindications

- Do not use over open wounds.
- Any known allergy and/or hypersensitivity to any material listed in product composition.

Warnings and Precautions

We recommend the initial fitting of this device be conducted by a suitably qualified healthcare professional who will advise the period of use.

- Carefully read all instructions and warnings prior to use.
- Follow all instructions to ensure proper performance of the device.
- Do not use if liniments, ointments, gels, creams or any other substances have been applied to the affected area.
- Do not re-use for another patient, doing so risks cross infection and can compromise product integrity.
- Perform regular skin and circulation checks, especially patients with diabetes, vascular deficiencies and neurological conditions.
- Should any adverse reactions occur, please cease use and contact your healthcare professional or provider of this product.
- The durability of the device may be compromised by certain factors, e.g., objects with sharp edges or damage to the hook and loop fasteners. To prevent this the hook and loop fasteners should always be fastened when the device is not being worn or when being washed.

Washing and Care Instructions

1. Close over all hook/loop straps.
2. Hand wash brace in cool water using mild detergent.
3. Rinse and dry flat.
4. Do not spin or tumble dry.
5. Do not use fabric conditioner at any time.
6. Check for any degradation that may affect use.
7. Periodically check hook/loop fasteners are operating securely.

Product Family Composition

Nylon, SBR, Polyester, Spandex.

Storage and Transport Conditions

Store in a cool, dry place out of direct sunlight, in the original packaging.

Recycling and Disposal

Packaging and constituent parts should be recycled or safely disposed in accordance with local or national laws.

Serious Incident

Report any serious incident to the manufacturer and the competent authority of the EU Member state, or country in which you reside.

MD MEDICAL DEVICE

1 SINGLE PATIENT - MULTIPLE USE

! CAUTION



Promedics Orthopaedics Ltd.
Port Glasgow Industrial Estate, Gareloch Road, Port Glasgow,
Scotland, UK, PA14 5XH
Tel: +44 (0) 1475 746 400 | Email: enquiries@promedics.co.uk
www.promedics.co.uk

CE **EC** **REP** **EMERGO EUROPE**
Westervoortsewijk 60,
6827 AT Arnhem
The Netherlands