



PRODUCT CATALOG



The
personal fit

IN ORTHOTIC FABRICATION
FOR PHYSICAL REHABILITATION

Dear Orfit customer,

A warm welcome to our product catalog, a place to learn about our orthotic fabrication materials and our company.

At Orfit Industries, we invent thermoplastic formulations and use innovative manufacturing techniques that offer the best solutions for you and your patients. Our 'light' products, our unique textile-like thermoplastics on a roll, and our all-time classics will help solve your most challenging patient cases.

However, not only performance counts. We also pay great attention to aesthetics by offering beautiful colors and a wide range of accessories to complete your orthoses and help increase patient adherence.

Finally, we provide training in orthotic techniques both at the level of inexperienced users and experienced clinicians.

Orfit is a family-owned company, with its owners close to the business. We are quick adapters to changing circumstances and navigate through the business with long-term implications in mind. We work hard not only to offer wonderful solutions to therapists and patients but also to ensure the well-being of our employees and stakeholders.

We hope you enjoy working with our products and invite you to turn to page 5 to discover our full product portfolio.

Yours sincerely,

A handwritten signature in blue ink, appearing to read 'Steven Cuypers', is written over a thin horizontal line.

Steven Cuypers
C.E.O, President

NEW

MULTIFIT® NS

ALL-ROUND
ORTHOTIC
FABRICATION
MATERIAL



NEW

ORFIGUARD

HIGH PRECISION
PROTECTION MASK



NEW

ISOFORCE


UNIQUE OUTRIGGER
SYSTEM FOR
STIFF PIP AND
IP JOINTS






NEW

STRAVA SPLINT BATH



| | |
|--|-----------|
| 1. WELCOME TO ORFIT INDUSTRIES | 3 |
| 2. ABOUT OUR PRODUCTS | 6 |
| 3. ORTHOTIC FABRICATION MATERIALS | 8 |
| Elastic Products | 8 |
| Orficast® Thermoplastic on a Roll | 8 |
| Orfit® Classic | 10 |
| Orfit® NS | 12 |
| Orfit® Colors NS | 14 |
| Orfit® Strips | 16 |
| Orfilight® NS | 18 |
| Aquafit™ NS (Orfit® Natural NS) | 20 |
| Tecnofit™ | 22 |
| Orfizip® Light NS | 24 |
| High Resistance to Stretch Products | 26 |
| Multifit® NS  | 26 |
| Orfit® Eco | 28 |
| Orfibrace™ NS | 30 |

| | |
|--|-----------|
| Precuts | 32 |
| Orfit® Classic precuts | 33 |
| Orfilight® NS precuts | 35 |
| Orfit® Eco precuts | 36 |
| Orfit® Colors NS precuts | 37 |
| Orfiguard  | 38 |
| Accessories | 40 |
| Allfit™ | 40 |
| Dynasyst™  Isoforce | 42 |
| Hook and Loop straps | 47 |
| Working tools | 48 |
| Padding and lining | 50 |
| Heating Devices | 51 |
| Water baths  | 51 |
| Dry heater | 53 |
| Heat gun | 54 |
| 4. PRODUCT OVERVIEW | 55 |
| 5. ACTIVITIES AND SERVICES | 74 |


ABOUT OUR PRODUCTS

These pages will introduce you to the different characteristics of our products.


Do you need help with selecting the right product for your use? Leave a message on orfit.com/contact

Perforations

No
Perforation



Macro
Perforation (1%)

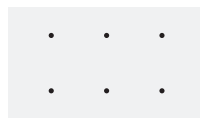

Maxi
Perforation (23.5%)




Mini
Perforation (3%)



Micro
Perforation (17%)



Micro Plus
Perforation (25%)

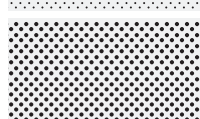

Perforated material allows for ventilation of the skin and makes the material lighter in weight.

Thicknesses

1.6 mm (1/16")



1.8 mm (1/14")



2.0 mm (1/12")



2.5 mm (3/32")



3.0 mm (1/8")



3.2 mm (1/8")



4.2 mm (3/16")



Colors

Neutral colors



White



Pearl White



Black



Dominant Black



Beige



Ivory

Metallic colors



Atomic Blue



Sonic Silver



Gold

Happy colors



Orange



Bright Pink



Dynamic Red



Blue



Ocean Blue



Hot Green



Violet

Non-stick properties

The majority of Orfit products carry a non-stick coating.

The non-stick (NS) coating allows for temporary bonding during fabrication. In addition, the NS coating protects the material from fingerprints and dirt.

Our coating is skin-friendly as it is water based and not based on solvents. The coating can be easily removed for applying outriggers.

Product classification: elasticity and resistance to stretch

Orfit products are classified by degree of elasticity and/or degree of resistance to stretch.

Resistance to stretch is the amount of resistance the material gives to being stretched when heated.

Materials with a high resistance to stretch require firm handling and are best suited for larger orthoses and for clients with high tone and/or contractures. Materials with low resistance to stretch require light handling to avoid over stretching.

Materials with high elasticity should be stretched over the body part for best conformity.

Precuts

Orfit precut orthoses are available in a variety of materials and designs. A precut orthosis will save you time and effort.

Simply activate in hot water, place over you patient's extremity and adjust only as needed.

ORFICAST® THERMOPLASTIC ON A ROLL

EXTREMELY COMFORTABLE FOR YOUR PATIENT

Available in: Orficast® 3 and 6 cm
Orficast® More 6, 12, 15 and 30 cm



Recommended applications

Perfect for traumatic and orthopedic conditions which require quick and easy orthoses.

Recommended as first aid product for quick immobilization and stabilization after injury.

Serial static orthoses to regain range of motion.

Recommended for children, new borns, clients with fragile skin and/or arthritic conditions.

Orficast is ideal for all types of finger orthoses and hand based wrist and thumb orthoses.

Orficast More is ideal for the fabrication of orthoses that go beyond the wrist: thumb-wrist opponens orthosis, De Quervain orthosis, radial/ulnar gutter orthosis and volar wrist cock-up orthosis.

Characteristics

Thermoplastic on a 3 meter roll

Comes in different widths, thicknesses and colors

Excellent elasticity in two directions

Adaptable when needed

Easy to cut

Self-adhesive capacity, especially when dry heated

Breathable material

How it works and looks

Material on a roll makes cutting a strip of material for orthotic fabrication easy and saves time. Easy to make orthoses even without a pattern.

A large variety in widths and thicknesses enables you to treat a wide range of clinical applications and to accommodate different orthotic designs. Orficast More allows, amongst others, for making a volar or dorsal slab.


Allows precise conformability to the patient's anatomy, for example in the web space.

Orthoses can be adapted when corrections are needed. Avoids waste of material.

Easy to cut and add small pieces to use as a reinforcement of the orthosis.

Dry heat is helpful for bonding corners and straps. Overlapping seams are not perceptible on the inside.

Excellent patient comfort and compliance: lightweight, allows ventilation and prevents maceration of the skin.

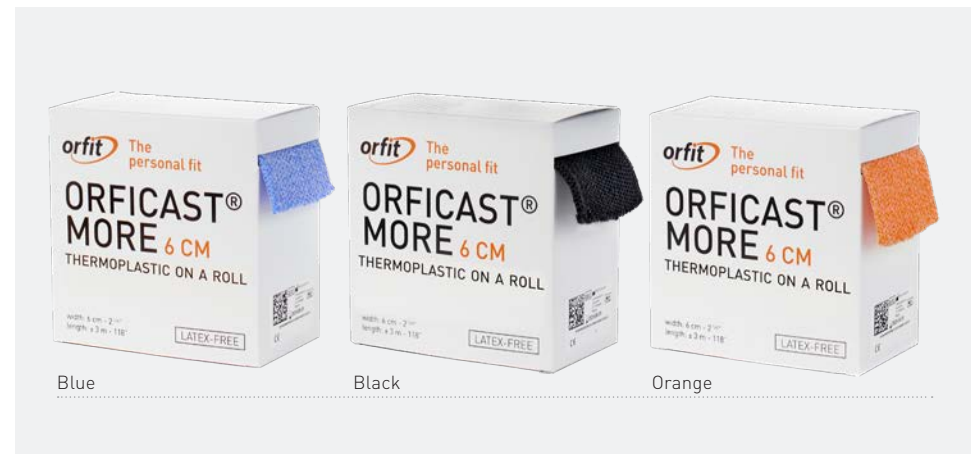

ORFIT® CLASSIC

HIGH PERFORMANCE HANDLING



Recommended applications

Unlimited potential for fabrication of very large to very small upper and lower extremity orthoses according to the thickness you choose.

“Soft” version:

Small finger, thumb, hand, wrist, wrist/thumb, elbow and toe orthoses and cervical collars.

“Stiff” version:

Large and long arm orthoses, AFO's and circumferentially applied fracture bracing.

Characteristics

Transparent when activated

Non-coated material

Excellent stretch (1000%) without tearing

Full elastic memory (100%)

Easy to trim edges

Available in high stretch ("soft") version
and low stretch ("stiff") version

How it works and looks

Transparency indicates when the material is ready to mold.
Allows for proper alignment and visual monitoring over anatomical landmarks.

Temporary adhesion during orthotic fabrication and reliable permanent bonding by using dry heat. Tackiness helps orthosis to stay in place during fabrication.
All attachments and outriggers can be applied with maximum ease.

Very precise molding that conforms to patient anatomy, for example in the hand arches and first web space.

Completely remoldable and easy to adjust after reheating locally. Avoids waste of material.

Extreme patient comfort.

"Soft" version has less resistance to stretch and more drape.
"Stiff" version has increased resistance to stretch and less drape.



ORFIT® NS

COATED FOR EASE OF USE



Recommended applications

Unlimited potential for fabrication of very large to very small orthoses.

Excellent material for the fabrication of hard shell pressure masks used to ensure excellent scar healing for the treatment of burned patients.

“Soft” version:

Small finger, thumb, hand, wrist, wrist/thumb, elbow and toe orthoses and cervical collars.

“Stiff” version:

Large and long arm orthoses, AFO's and circumferentially applied fracture bracing.

Characteristics

Transparent when activated

Coated, non-stick material

Excellent stretch (1000%) without tearing

Full elastic memory (100%)

Easy to trim edges

Available in high stretch ("soft") version
and low stretch ("stiff") version

How it works and looks

Transparency indicates when the material is ready to mold.
Allows for proper alignment and visual monitoring over anatomical landmarks.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication.
Reliable permanent bonding is possible after removing the coating and using dry heat.
All attachments and outriggers can be applied with maximum ease.

Very precise molding that conforms to patient anatomy, for example in the hand arches and first web space. Accurate molding improves correct dispersion of skin pressure.

Completely remoldable and easy to adjust after reheating locally.
Avoids waste of material.

Extreme patient comfort.

"Soft" version has less resistance to stretch and more drape.
"Stiff" version has increased resistance to stretch and less drape.

ORFIT® COLORS NS

EXCITING COLORS FOR YOUR PATIENTS



Recommended applications

Unlimited potential for fabrication of very small to very large orthoses according to the thickness you choose. Small finger, thumb, hand, wrist, wrist/thumb, elbow and toe orthoses, cervical collars and circumferentially applied fracture braces.

Small strips can be used as nice and colorful edging strip.

Characteristics

Unique painter's palette of metallic and bright colors with a nice high gloss and shiny surface

Coated, non-stick material

Excellent stretch (1000%) without tearing

Excellent elastic memory (100%)

Easy to trim edges

How it works and looks

Wide range of colors to meet every patient's personal preferences, certainly children and sport athletes. Combined with a lot of other beneficial factors, this will improve patient compliance.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Very precise molding that conforms to patient anatomy, for example in the hand arches and first web space. Accurate molding improves correct dispersion of skin pressure.

Completely remoldable and easy to adjust after reheating locally. Avoids waste of material.

Extreme patient comfort.



ORFIT® STRIPS

QUICK AND EASY TO USE



Recommended applications

Quick and easy molding of small corrective finger and hand orthoses, such as anti-swan neck orthosis and anti-claw hand orthosis.

Can be used to fabricate handles and larger grips for everyday objects to increase functional independence in individuals with limited finger motion.

Adheres to other orthoses to increase rigidity.

Characteristics

Precut strip of Orfit Colors NS

Coated, non-stick material

Smooth edges cut by machine

Available in two widths and two thicknesses

How it works and looks

Less waste of material. Lightweight and beautiful orthoses.
Choice of color according to patient preferences.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Quick finish of the orthosis.

Choice of material thickness according to the required resistance to force.



ORFILIGHT® NS

UP TO 30% LIGHTER THAN OTHER MATERIALS

Available in: Orfilight® NS, Orfilight® Black NS
and Orfilight® Atomic Blue NS



Recommended applications

Excellent material choice for all types of small, medium and large orthoses.

Ideal material to fabricate post-surgical orthoses.

Especially recommended for patients with loss of strength and/or fragile and sensitive skin, for example, the elderly, pediatric patients or for arthritic conditions including rheumatoid arthritis and others.

Lightweight base for all outrigger attachments.

Characteristics

Weighs 25% to 30% less than other materials

Coated, non-stick material

Excellent stretch, moderate elasticity

Maximum stretch without tearing

Moderate memory

Soft velvet surface touch and easy to trim edges

How it works and looks

Allows to make very lightweight and comfortable orthoses.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Easy molding for orthosis to anatomy.

Easy to handle and perfect control of material during molding. Material can be stretched sufficiently in order to position the material correctly and hold it in place during orthotic fabrication.

Completely remoldable and easy to adjust after reheating locally. Avoids waste of material.

Patient comfort leading to increased patient compliance with wearing schedule.



AQUAFIT™ NS (ORFIT® NATURAL NS)

THE BEAUTY OF CRYSTAL WHITE



Recommended applications

Unlimited potential for fabrication of very large to very small upper and lower extremity orthoses according to the thickness you choose.

Excellent material for the fabrication of hard shell pressure masks used to ensure excellent scar healing for the treatment of burned patients.

“Soft” version:

Small finger, thumb, hand, wrist, wrist/thumb, elbow and toe orthoses and cervical collars.

“Stiff” version:

Large and long arm orthoses, AFO's and circumferentially applied fracture bracing.

Characteristics

Transparent when activated

Coated, non-stick material

Excellent stretch (1000%) without tearing

Full elastic memory (100%)

Easy to trim edges

Available in high stretch ("soft") version and low stretch ("stiff") version

How it works and looks

Transparency indicates when the material is ready to mold.
Allows for proper alignment and visual monitoring over anatomical landmarks.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication.
Reliable permanent bonding is possible after removing the coating and using dry heat.
All attachments and outriggers can be applied with maximum ease.

Very precise molding that conforms to patient anatomy, for example in the hand arches and first web space. Accurate molding improves correct dispersion of skin pressure.

Completely remoldable and easy to adjust after reheating locally. Avoids waste of material.

Extreme patient comfort.

"Soft" version has less resistance to stretch and more drape.
"Stiff" version has increased resistance to stretch and less drape.



TECNOFIT™

STRONG AND RIGID SUPPORT



Recommended applications

Unlimited potential for the fabrication of all size orthoses according to the thickness you choose with extra care and prevention of pressure points.

Thumb, hand, wrist, wrist/thumb and elbow orthoses, cervical collars and body jackets.

Large and long arm orthoses, AFO's and fracture bracing.

Excellent as a base for different attachments.

Characteristics

Transparent when activated

Coated, non-stick material

Excellent stretch (300%) without tearing

Excellent elastic memory

Easy to trim edges

How it works and looks

Transparency indicates when the material is ready to mold.
Allows for proper alignment and visual monitoring over anatomical landmarks.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Very precise molding that conforms to patient anatomy, for example in the hand arches and first web space. Accurate molding improves correct dispersion of skin pressure.

Completely remoldable and easy to adjust after reheating locally. Avoids waste of material.

Extreme patient comfort.

ORFIZIP® LIGHT NS

ACTIVATE, MOLD AND ZIP



Recommended applications

Lightweight circumferential orthoses for immobilization and protection following trauma and/or surgery.

Foamy surface is especially beneficial for patients with sensitive and/or fragile skin.

Use with mobilization orthoses: outriggers are easily bonded to the surface and the base orthosis will not migrate distally.

Offers a snug fit without straps.

Characteristics

Precut wrist orthosis with double-sided zipper

Excellent stretch

Moderate elastic memory

Coated, non-stick material

Can accommodate all wrist sizes

Quick and easy one-step application

Increased patient comfort and compliance

How it works and looks

Each precut can be used to fabricate either a left or right hand orthosis.

Very precise molding conforms to patient anatomy, for example in the hand arches and first web space. Accurate molding improves correct dispersion of skin pressure.

Can be reactivated and remolded.
Easy to adjust after reheating locally, which avoids waste of material.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Available in two thicknesses and three sizes (S, M, L) to accommodate all patients.

Time-saving (no additional attachments needed).
Length of zipper can be adjusted (zipper must be trimmed and properly sealed).

Washable, extra lightweight material with a soft velvet surface touch and optimal skin ventilation. Easy to put on and take off.



MULTIFIT® NS

EMBRACE DIVERSITY

Available in: Pearl White and Black



Recommended applications

Multifit NS is user friendly for the fabrication of all types of orthoses.

Use Multifit NS 3.2 mm for large orthoses and patients that require extra support: resting hand orthosis, volar wrist orthosis, elbow orthosis, use as a base for outrigger components and contracture management.

Characteristics

High rigidity

High resistance to stretch and moderate drape

Minimal memory

Coated, non-stick material

How it works and looks

Ideal material for the fabrication of rigid and supportive orthoses. It is an excellent base material for all outrigger and resting hand orthoses.

The perfect balance between drape and resistance to stretch. The material holds its shape perfectly against increased tone, and it readily contours to the fingers, hands and wrist using mild to moderate pressure.

The fabricated orthosis will hold its shape when heated for minor adjustments. Small corrections are easy to make after spot-heating the orthosis with a heat gun.

The non-stick coating prevents accidental bonding and allows for molding directly over bandages and cotton linings in post-surgical conditions. The coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.



ORFIT® ECO

THE PERFECT SOLUTION FOR TIGHT BUDGETS

Available in: Orfit® Eco and
Orfit® Eco Black NS



Recommended applications

Excellent rigidity for large orthoses requiring full support:

- > long arm and elbow orthoses
- > knee and foot orthoses
- > body jackets
- > back supports

Specifically recommended to fabricate anti-spastic orthoses.

Characteristics

Orfit Eco is a tacky, non-coated material

Orfit Eco Black NS is a non-stick, coated material (NS version)

Offers excellent rigidity and excellent resistance to fingerprints

High resistance to stretch and low drape

Low to moderate elastic memory

High quality recycled material

Beautiful smooth surface finish

How it works and looks

Temporary adhesion during orthotic fabrication and reliable permanent bonding by using dry heat. Tackiness helps orthosis to stay in place during fabrication. All attachments and outriggers can be adhered easily.

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Easy handling and control of the material during molding. Allows to apply pressure during the fabrication of the orthosis.

Very strong material which holds its shape perfectly against increased tone. The combination of the molding capacity and a precise pattern design allows to mold and conform to patient anatomy.

Orthoses are completely remoldable. Avoids waste of material.

Ecological material.

Excellent patient comfort.



ORFIBRACE™ NS

BACK TO ACTION, AS SOON AS POSSIBLE



Recommended applications

All large and rigid orthoses.

Trunk corsets, ankle supports, elbow, forearm and wrist orthoses.

Permanent bracing.

Ideal to mold onto positive plaster molds protected by a cotton lining.

Characteristics

Coated, non-stick material

Excellent rigidity and excellent resistance to fingerprints

High resistance to stretch and low drape

Moderate memory

Easy to trim and smooth edges

Smooth surface finish

How it works and looks

No accidental bonding. Allows for molding directly over bandages and cotton linings in post-surgical conditions. Coating allows for temporary bonding during fabrication. Reliable permanent bonding is possible after removing the coating and using dry heat. All attachments and outriggers can be applied with maximum ease.

Allows to apply pressure during the fabrication of the orthosis. Excellent for all mechanical attachments.

Easy handling and control of material during molding.



Remoldable. Orthosis can be easily adjusted after local reheating.

Extreme patient comfort.

Comfortable and pleasant surface feel, which may contribute to patient compliance.



ORFIT PRECUTS

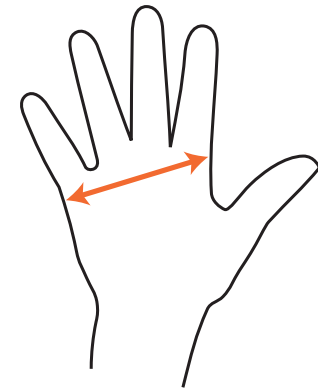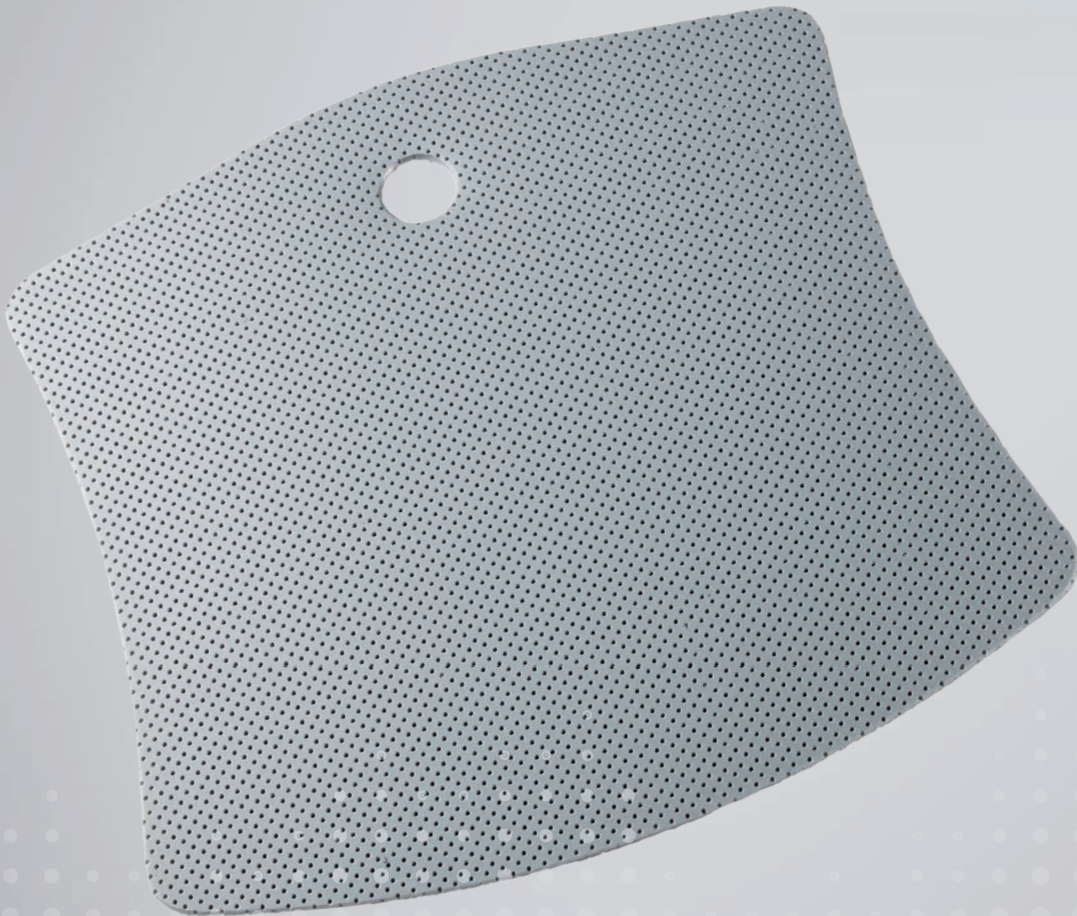

SIZING

Measure the total hand width at the MCP joints of the 2nd to 5th finger.

60-70 mm (2 ¹/₅" - 2 ³/₄") = **SMALL**

80-90 mm (3 ¹/₈" - 3 ¹/₂") = **MEDIUM**

100-110 mm (4" - 4 ³/₁₀") = **LARGE**

ORFIT® CLASSIC PRECUTS

**GAUNTLET THUMB POST ORTHOSIS**

Orfit® Classic


1.6 mm (1/16") micro perfo, soft

This orthosis gives an extremely comfortable immobilization and stabilization of the complete thumb column. It can serve for the treatment of De Quervain tenosynovitis, rheumatoid arthritis, basal joint osteoarthritis, thumb sprain and pre- or post-operative thumb management.

35810 (small)

35811 (medium)

35812 (large)

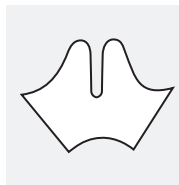**COMPLETE THUMB ORTHOSIS**

Orfit® Classic


2.0 mm (1/12") micro perfo, soft

Intended for acute cases. The dorsal part at the level of the thumb prevents hyperextension of the distal phalanx. This part can easily be cut away when necessary.

35900 (small)

35901 (medium)

35902 (large)

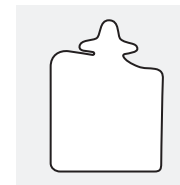**WRIST + THUMB ORTHOSIS**

Orfit® Classic


2.0 mm (1/12") micro perfo, soft

A traditional wrist and thumb immobilization orthosis for very acute joint and tendon inflammations. This orthosis gives full comfort thanks to its circumferential design and long lever.

35814 (small)

35815 (medium)

35816 (large)

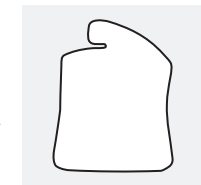**GAUNTLET IMMOBILIZATION ORTHOSIS**

Orfit® Classic


2.0 mm (1/12") micro perfo, soft

A complete immobilization orthosis for the wrist. Due to its lightweight and small size, it is very comfortable, very supportive and also extremely strong and durable.

35830 (small)

35831 (medium)

35832 (large)

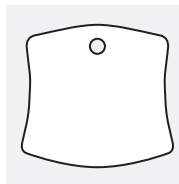**DORSAL WRIST COCK-UP ORTHOSIS**

Orfit® Classic

3.2 mm - 1/8" non perfo, stiff

An excellent alternative to volar wrist cock-up orthosis as this design allows for increased function while maintaining the wrist in extension. This orthosis design is suitable as a base for the attachment of low profile outriggers. Available in four sizes for all hand sizes from children to adults.

35818 (x-small)

35820 (small)

35821 (medium)

35822 (large)

**INTRINSIC RESTING HAND ORTHOSIS**

Orfit® Classic


3.2 mm - 1/8" non perfo

An orthotic design for management of a variety of hand problems.

35850 (small)

35851 (medium)

35852 (large)

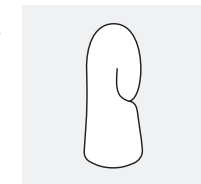

ORFIT® CLASSIC PRECUTS



ANTI-SPASTIC ORTHOSIS + THUMB PIECE

Orfit® Classic


3.2 mm - 1/8" non perfo, stiff

This orthosis design offers an easy way to manage the patient with spasticity of the upper extremity. The volar-dorsal design and separate thumb piece make orthotic fabrication easier. The orthosis maintains the wrist and fingers in extension.

35870 (small)

35871 (medium)

35872 (large)

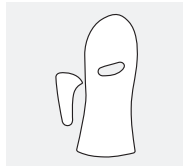

RHEUMATOID ARTHRITIS RESTING ORTHOSIS

Orfit® Classic


3.2 mm - 1/8" mini perfo

A specific design for preventive and corrective orthotic fabrication of the stereotype rheumatoid hand. Both antero-posterior and latero-lateral deformities of wrist and fingers, as well as thumb deformities, are easily realigned.

35840 (small)

35841 (medium)

35842 (large)

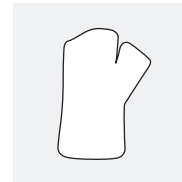

ULNAR DRIFT ORTHOSIS


Orfit® Classic

3.2 mm - 1/8" non perfo, stiff

This unique pattern offers the most realistic and functional ulnar deviation orthosis for patients with early stages of rheumatoid arthritis and ulnar drift of the digits.

35817KL (left)

35817KR (right)

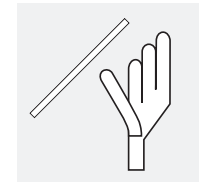

ORFILIGHT® NS PRECUTS

**GAUNTLET THUMB POST ORTHOSIS**

Orfilight® NS


1.6 mm (1/16") micro perfo

This orthosis gives an extremely comfortable immobilization and stabilization of the complete thumb column. It can serve for the treatment of De Quervain tenosynovitis, rheumatoid arthritis, basal joint osteoarthritis, thumb sprain and pre- or post-operative thumb management.

35810/L (small)

35811/L (medium)

35812/L (large)

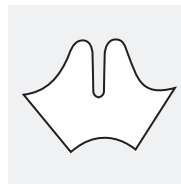**DORSAL WRIST COCK-UP ORTHOSIS**

Orfilight® NS

3.2 mm - 1/8" non perfo

An excellent alternative to volar wrist cock-up orthosis as this design allows for increased function while maintaining the wrist in extension. This orthosis design is suitable as a base or the attachment of low profile outriggers.

35820/L (small)

35821/L (medium)

35822/L (large)

**GAUNTLET THUMB POST ORTHOSIS**

Orfilight® Black NS


1.6 mm (1/16") micro perfo

This orthosis gives an extremely comfortable immobilization and stabilization of the complete thumb column. It can serve for the treatment of De Quervain tenosynovitis, rheumatoid arthritis, basal joint osteoarthritis, thumb sprain and pre- or post-operative thumb management.

35810Z/L/NS (small)

35811Z/L/NS (medium)

35812Z/L/NS (large)

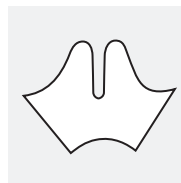**WRIST + THUMB ORTHOSIS**

Orfilight® Black NS

2.5 mm (3/32") micro perfo

A traditional wrist and thumb immobilization orthosis for very acute joint and tendon inflammations. This orthosis gives full comfort thanks to its circumferential design and long lever.

35814Z/L/NS (small)

35815Z/L/NS (medium)

35816Z/L/NS (large)



ORFIT® ECO PRECUTS



DORSAL WRIST COCK-UP ORTHOSIS

Orfit® Eco

3.2 mm (1/8") non perfo

A most functional orthosis for most wrist complaints, especially for flaccid palsy. Very suitable as a basic orthosis to receive low profile extension outriggers.

35820/ECO (small)

35821/ECO (medium)

35822/ECO (large)



ANTI-SPASTIC ORTHOSIS + THUMB PIECE

Orfit® Eco Black NS

3.2 mm (1/8") mini perfo

The purpose of this orthosis is to maintain the long flexor muscles stretched in order to inhibit the hypertonicity.

35870Z/ECO/NS (small)

35871Z/ECO/NS (medium)

35872Z/ECO/NS (large)



INTRINSIC ANTI-SPASTIC HAND ORTHOSIS

Orfit® Eco Black NS

3.2 mm (1/8") mini perfo

This orthosis has a radial bar for added strength on the radial side of the orthosis, ideal for patients with spasticity.

35910Z/ECO/NS (small)

35911Z/ECO/NS (medium)

35912Z/ECO/NS (large)



ORFIT® COLORS NS PRECUTS

**GAUNTLET THUMB POST ORTHOSIS**

Orfit® Colors NS


2.0 mm [1/12"] micro perfo

This orthosis gives an extremely comfortable immobilization and stabilization of the complete thumb column. It can serve for the treatment of De Quervain tenosynovitis, rheumatoid arthritis, basal joint steoarthritis, thumb sprain and pre- or post-operative thumb management.

35810MB/NS (small)

35811MB/NS (medium)

35812MB/NS (large)

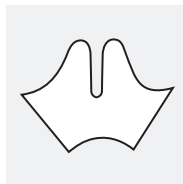**WRIST + THUMB ORTHOSIS**

Orfit® Colors NS

2.0 mm [1/12"] micro perfo

A traditional wrist and thumb immobilization orthosis for very acute joint and tendon inflammations. This orthosis gives full comfort thanks to its circumferential design and long lever.

35814Z/NS (small)

35815Z/NS (medium)

35816Z/NS (large)

**GAUNTLET IMMOBILIZATION ORTHOSIS**

Orfit® Colors NS


2.0 mm [1/12"] micro perfo

A complete immobilization orthosis for the wrist. Due to its lightweight and small size, it is very comfortable, very supportive and also extremely strong and durable.

35830MG/NS (small)

35831MG/NS (medium)

35832MG/NS (large)

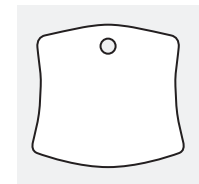

Orfit Colors NS precuts are available in

- > 3 designs
- > 4 colors (Atomic Blue, Ocean Blue, Sonic Silver and Dominant Black)
- > 3 sizes

See Product Overview for the full list of article numbers.

Orfit precuts do not include accessories, such as hook and loop straps or edging strips.

ORFIGUARD

HIGH PRECISION PROTECTION MASK



Product Characteristics

Orfiguard comes as a complete set with strapping material and the mask. With the one size fits all thermoplastic precut mask, no specific measurements need to be taken for sizing.

The material used for this mask is Orfit Colors NS in 2.0 mm, a durable, strong and long lasting material. The mask is completely remoldable and easy to adjust after reheating locally.

Orfiguard can be trimmed for younger patients.



Orfiguard shields the face after a facial injury, fracture or other condition which requires protection of the nose, nose bridge and/or eye sockets.

ALLFIT™

ALWAYS FITS WHATEVER YOUR PATIENT NEEDS



Recommended applications

Change the shape of any appliance or utensil handle.

Add additional components to orthoses.

Use as first web spacer or finger separator.

Improve patient grip by enlarging, enhancing or reshaping the objects one needs to hold and use during daily life activities, for example, knives, forks and spoons.

Characteristics

Small pellets of thermoplastic material

Convenient sized bags for multiple applications

Quick activation time (1 minute)

Transparent when activated, white when cooled

Pellets simply stick together uniformly and stick very well to other surface objects.

How it works and looks

No need to cut small scraps from larger sheets.

Easy to preserve.

Fast and easy to use and no need for special training.

Indicates when the material is ready to mold. Pleasing neutral color.

Short melting time and quick adherence around objects.

Easy to mold into a desired shape. Smooth surface all around for increased patient comfort.

Simply pour the pellets into a small container of hot water 65°C (149°F). When they have transformed into a transparent mass, Allfit is ready to use.



DYNASYST™

MOBILIZATION ORTHOSES WITH ORFITUBE, COILS AND SPRINGS, AND ISOFORCE



Recommended applications

Orfitube

- Mobilization orthoses, for example post-operative orthoses after tendon surgery or to recover range of motion
- Fabrication of postural orthoses

Coils and Springs

- Dynamic orthoses for nerve palsies
- Dynamic PIP joint extension orthoses for extensor lag at the PIP joint

Isoforce

- Lack of full passive PIP joint extension
- PIP joint flexion contracture
- Dynamic and static progressive orthoses

Characteristics

Lightweight outrigger attachments

Patient pleasing aesthetic appearance

Low profile design

Large choice between Orfitubes and joint spring coils

Fast and easy fabrication for outrigger attachments

Outrigger system specifically designed to address the problem of the stiff PIP joint

How it works and looks

Increased patient comfort and wearing schedule.

May lead to increased patient compliance with protocol.

Spring coils are placed at the lateral sides of the joints.
Crossbar systems create a perpendicular traction force on the involved segment.


Orfitube Adjusters allow for easy and correct alignment of finger loops and easy to adjust angles.
Simplified adjustment and re-adjustment of correct angles of alignment and applied forces.

Time-saving.

Isoforce ensures a constant force transmission and a constant force application angle throughout the complete PIP joint range of motion. Isoforce maintains a 90 degrees force application angle regardless of changes in joint position.



ORFITUBE

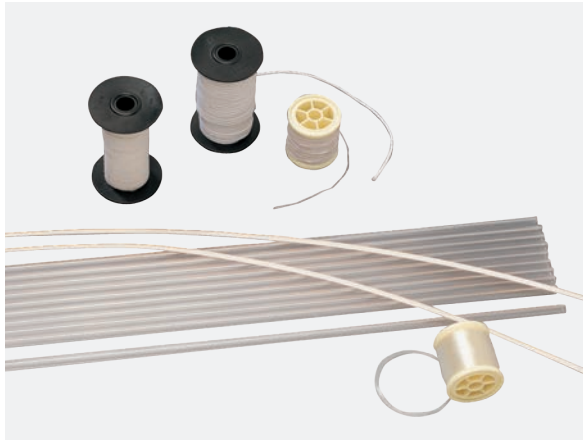

Art. N° 35311 - ORFITUBE™ KIT

Contains:

- > 10 tubes of 330 mm (13") long
- > 3 different types of elastic thread with different degrees of elastic strength (on a spool).
- > 1 spool polyester monofilament for an optimal gliding effect (if necessary).
- > 2 bending wires.




Art. N° 35305 - ORFITUBE™ ADJUSTER KIT

The Orfitube Adjuster puts you in full control of longitudinal and lateral position of Orfitubes. But there is more than one angle from which you can look at this product. The compact design allows you to fix one individual adjuster for each finger. The adjusters are light-weight, very cost-effective and more elegant than ever before.



Art. N° 35318 - ORFITUBE™ ENDS

20 pcs.
Chromium-plated tube ends, fit on all Orfitubes, reducing friction of the elastic thread with the tube edge to a minimum.

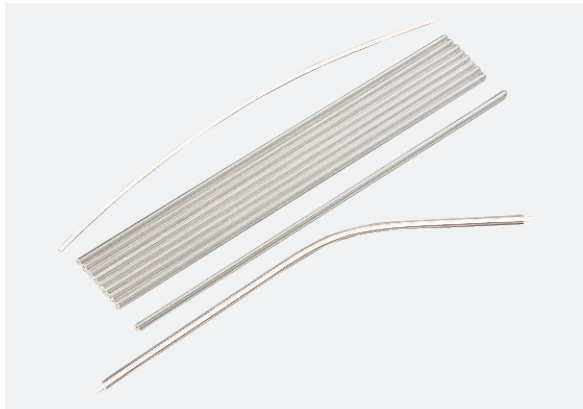

Art. N° 35310 - ORFITUBES™ + BENDING WIRES


10 tubes of 330 mm (13") long, incl. 2 bending wires.



Art. N° 35304 - SINGLE ORFITUBE™ ADJUSTER - 5 PCS. +

ALLEN WRENCH


5 single Orfitube adjusters + 1 Allen wrench

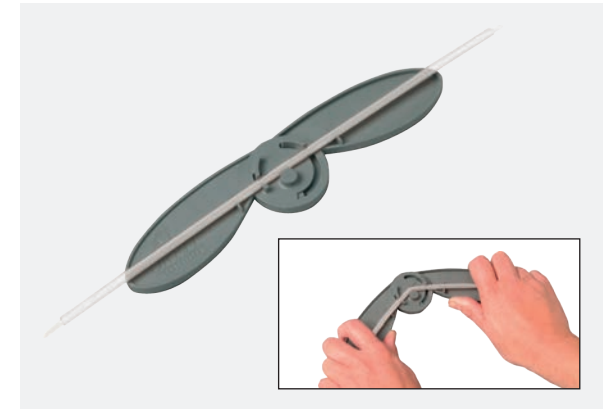

Art. N° 35308 - ORFITUBE™ BENDING TOOL


For easy bending of the strong Orfitubes without stress on the thumbs.

COILS AND SPRINGS

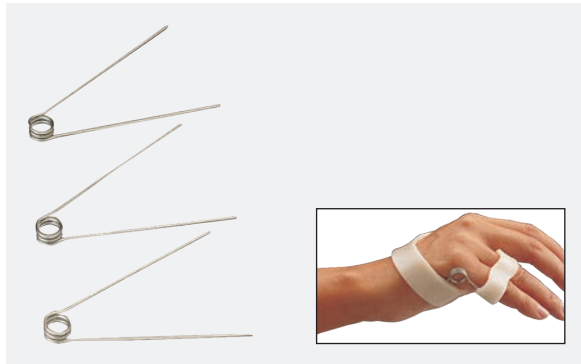**Art. N° 35801 - KNUCKLE BENDER COIL SPRING**

10 pcs. (Gauge 1 mm [1/24"] - coil diameter 10 mm [5/12"])

Designed to assist the intrinsic muscle weakness, e.g. with ulnar nerve palsy. The spring can also be turned upside down to increase spring power. Torque force at full finger extension: ca. 1000 g.cm.

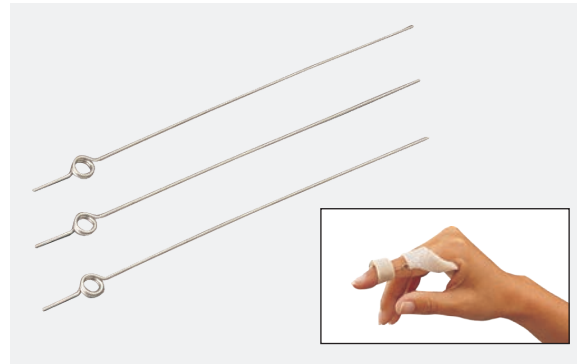**Art. N° 35800 - FINGER COIL SPRING**


10 pcs. (Gauge 0.7 mm [1/37"] - coil diameter 5 mm [1/5"])

The short leg is bent and anchored in the dorsal bridge of the orthosis. The long leg is bent and cut to length. The spring power is adjusted by altering the coil design.

**Art. N° 35802 - FINGER EXTENSION ASSIST**

10 pcs. (Gauge 1 mm [1/24"] - coil diameter 4 mm [1/6"])

The short leg is secured between two layers of Orfit Classic. The long leg is bent to size. At total finger flexion the spring is tightened up to a torque force of ca. 2000 g.cm. The spring is especially designed for finger extension assists in combination with the Orfit orthosis Art. N° 35822/ECO.

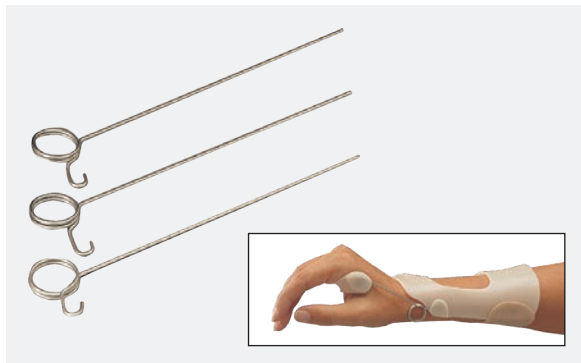**Art. N° 35804 - WRIST EXTENSION COIL SPRING**

10 pcs. (5 left + 5 right)

(Gauge 1.5 mm [1/16"] - coil diameter 15 mm [3/5"])

Two types supplied: coiled (left or right)

Torque force when tightened over an angle of 130°: 5000 g.cm.

**RADIAL WRIST EXTENSION SPLINT**

Art. N° 35865K (small)

Art. N° 35866K (medium)

Radial nerve patients need two strong assisting extension coils.

When recovering, the coils help to rebalance muscle power.

Set for fabrication of one orthosis includes:

> 2 Orfit Classic precuts (3.2 mm [1/8"] non perfo stiff).

> 2 lateral coil springs: gauge diameter 1.5 mm [1/16"]

Coil diameter 15 mm [3/5"] - Torque force: 5000 g. cm.

**Art. N° 35805 - SPRING WIRE**

Length: 15 m [591"]

Thickness: 1 mm [1/24"]

For making your own coil springs and outriggers.

ISOFORCE



ISOFORCE

Art. N° 35322 (small)

Art. N° 35323 (medium)

Art. N° 35324 (large)

A unique outrigger system specifically designed to address stiff proximal interphalangeal (PIP) and interphalangeal (IP) joints. It can be used for the creation of dynamic and static progressive orthoses.

THREADS



Art. N° 35315 - POLYESTER MONOFILAMENT THREAD

Spool of 75 m [2953"]

Used to increase the gliding of the elastic thread in the Orfitube, e.g. when the Orfitube is locally kinked and when the length of the elastic is not determined.



ELASTIC THREAD

Spool of 10 m [394"]

Stretchability: at least 200%. No ageing after 6 weeks.

Maximum stretching.

Art. N° 35312: Max. resistance \pm 200 g.

For supporting weakened finger muscles.

Art. N° 35313: Max. resistance \pm 400 g.

For most dynamic finger orthoses.

Art. N° 35314: Max. resistance \pm 600 g.

For dynamic wrist extension orthoses or for very persistent finger contractures.



HOOK AND LOOP STRAP 25 MM (1")

Available in a standard and self-adhesive version in different colors.




HOOK AND LOOP STRAP 50 MM (2")

Available in a standard and self-adhesive version in different colors.



DOUBLE SIDED HOOK OR LOOP STRAP

25 mm (1") width x 25 m (82') length

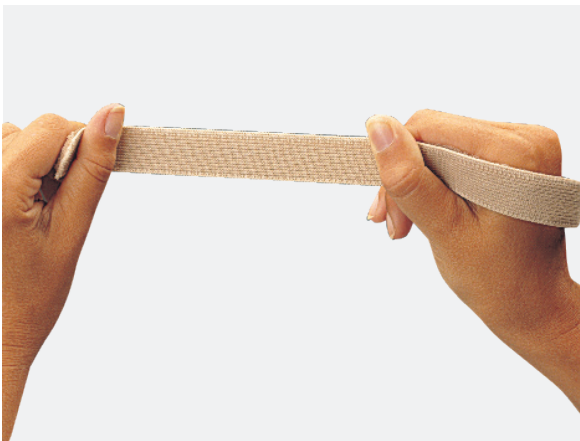

ELASTIC LOOP TAPE

Available in different widths and colors.



SPLINT STRAP

Available in different lengths.

See full product range
on the product overview pages.



Art. N° 35061 - SHARP SCISSORS - CURVED (14 CM)

Slim design excellent for trimming orthoses. The stainless steel blades are curved to make cutting around curves easy.

Surgical grade instrument.

Scissor length: 14 cm - 5 33/64"



Art. N° 35065 - STRONG SCISSORS

Strong cogged pair of scissors with return spring to cut sheets from 2 mm (1/12") to 4.2 mm (3/16").



Art. N° 35066 - GINGER SHEARS

Voted best scissors in the world by the majority of the therapists.

Scissor length: 20 cm - 7 7/8"



Art. N° 35067 - BANDAGE SCISSORS

Ideal for removing and trimming of thermoplastic splinting material.

Stainless steel blades with blunt tips to provide safety during removal.



Art. N° 35140 - SPRING COIL PLIERS

Full round nosed pliers for turning spring wire into coil with a diameter from 5 mm - 1/5" to 10 mm - 2/5".



Art. N° 35134 - CHROME PLATED EYELET SETTER

To set eyelets in leather, fixation straps and Orfit 1.6 mm - 1/16" and 2 mm - 1/12".

Art. N° 35137: 200 eyelets of 4 mm (1/6" diameter)



Art. N° 35307 - DEBURRING KNIFE

To smoothen orthosis edges. Ideal to finish orthoses made of non-perforated material in a fast and efficient way. Two blades are inside the handle.

Art. N° 35307B10: 10 spare blades with straight cut

Art. N° 35307B20: 10 spare blades with oblique cut



Art. N° 35105 - DYNAMOMETER

Accurate dynamometer, max. 600 g.

Ideal to measure static and dynamic orthosis tractions.

A ring is attached to one end, a clip to the other.



Art. N° 35003G - YELLOW MARKER

This special marker allows drawing a pattern on the splinting material that is indelible even in water, but is easily wiped off.

Art. N° 35003G/1 - Refill Yellow Marker



Art. N° 35003Z - BLACK MARKER

This special marker allows drawing a pattern on the splinting material that is indelible even in water, but is easily wiped off.

Art. N° 35003Z/1 - Refill Black Marker



Art. N° 35003W - WHITE MARKER

This special marker allows drawing a pattern on the splinting material that is indelible even in water, but is easily wiped off.

Art. N° 35003W/1 - Refill White Marker



Art. N° 35110 - LUXOFOAM NON PERFORATED

Luxofoam is a self-adhesive soft PVC foam used for padding of small surfaces. It does not detach when heated in water, and can therefore be used for prepadding orthosis applications.

Sizes: 49 cm x 97 cm x 3 mm (19 3/10" x 38 1/5" x 7/64")



Art. N° 35112 - LUXOFOAM PERFORATED

Sizes: 49 cm x 97 cm x 3 mm (19 3/10" x 38 1/5" cm x 7/64")



Art. N° 35111 - ORFILASTIC

Orfilastic is a terrycloth lining material with a backing of self-adhesive open cell foam. Used for covering large orthosis surfaces to prevent and absorb perspiration.

Easily removable with lukewarm water.

Sizes: 50 cm x 1 m x 3 mm (20" x 39 2/5" x 7/64")



Art. N° 35117 - ORFILASTIC BLACK

Orfilastic Black is a terrycloth lining material with a backing of self-adhesive open cell foam. Used for covering large orthosis surfaces to prevent and absorb perspiration. Easily removable with lukewarm water.

Sizes: 50 cm x 1 m x 3 mm (20" x 39 2/5" x 7/64")



Art. N° 360201N.ST1 - EDGING STRIP NATURAL COLOR

Thin low temperature thermoplastic lining material.

Sizes: 0,8 mm x 20 mm x 3000 mm (0 1/32" x 3/4" x 118 1/10")



Art. N° 35116 - SOFT VELVET EDGING STRIP BLACK

Thin, soft, stretchable, elastic and self-adhesive textile lining material.

Sizes: 1 mm x 25 mm x 2500 mm (1/25" x 1" x 8' 2.5")

NEW

STRAVA SPLINT BATH

**Art. N° 35123/230EU – STRAVA SPLINT BATH**

Inner dimensions: 480 L x 300 W x 150 H mm

Temperature range: 60°C - 82°C

Water content: approximately 20 liter



Automatic, high-heat sanitation cleaning

Continuous UV-C for effective sterilization

Only 25 minutes to reach set water temperature

Uniform heat distribution

Automatic on/off

Touch display with custom preset heating temperatures and times

Visual timer with audible alert

WATER BATHS



Art. N° 32502/230EU - WATER BATH

Water bath in premium stainless steel.

Inner dimensions: 500 L x 300 W x 103 H mm

Temperature range: 60°C - 90°C

Water content: approximately 10 liter



Art. N° 35097/115/1400 (115V) - LARGE WATER BATH (SUSPAN 1)

Art. N° 35097/230/1400 (220V) - LARGE WATER BATH (SUSPAN 1)

Large water bath in stainless steel.

Inner dimensions: 740 L x 540 W x 100 H mm

Temperature range: 60°C - 90°C

Water content: approximately 30 liter



Integrated faucet facilitates regular draining and replacing of water.



Excellent temperature control and temperature uniformity with a digital display and timer.



Grid with high and accessible handles provides easy and safe loading and removal of products.

DRY HEATER



Heat surface dimensions: 438 L x 318 W mm

Temperature range: 60°C - 85°C

Delivered with detachable power cable, teflon sheet, heating guide and user manual.

Create a convenient workplace with the Orfit Dry Heater and enjoy the pleasure of a water-free and fast working method.



Full power at your fingertips
Fast activation of your orthotic fabrication material.



Smart casing
Portable case that makes your work easy on the go.



Effective heating control
Excellent temperature control and temperature uniformity with a digital display and integrated timer.

DRY HEATER

| Art. N° | Voltage | Region/Country |
|-------------|---------|----------------|
| 32503/230EU | 230V | EU |
| 32503/230UK | 230V | UK |
| 32503/230CH | 230V | CH |

HEAT GUN

Art. N° 35353



Electrical heat gun.

Temperature range: variable heat settings
between 50°C - 650°C.

Adjustable heat flow: 250 - 500 L/min.

230V - 2000 Watt.

**Ideal to heat precise areas
on small and large thermoplastic
orthoses to attach outriggers
and straps or to make detailed
corrections.**



PRODUCT OVERVIEW

Consult the next pages for an overview of all products, including the article numbers.

PRODUCT OVERVIEW



ORFICAST® 3 CM

ORFICAST® 3 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|-------|---------|
| 1 roll | 3 cm | Blue | 4030 |
| | | Black | 4030Z |
| 12 rolls | 3 cm | Blue | 4031 |
| | | Black | 4031Z |



ORFICAST® 6 CM

ORFICAST® 6 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|--------|---------|
| 1 roll | 6 cm | Blue | 4032 |
| | | Black | 4032Z |
| | | Orange | 4032OR |
| 6 rolls | 6 cm | Blue | 4033 |
| | | Black | 4033Z |
| | | Orange | 4033OR |



ORFICAST® MORE 6 CM

ORFICAST® MORE 6 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|--------|---------|
| 1 roll | 6 cm | Blue | 4035 |
| | | Black | 4035Z |
| 6 rolls | 6 cm | Orange | 4035OR |
| | | Blue | 4036 |
| | | Black | 4036Z |
| | | Orange | 4036OR |



ORFICAST® MORE 12 CM

ORFICAST® MORE 12 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|-------|---------|
| 1 roll | 12 cm | Blue | 4037 |
| | | Black | 4037Z |
| 6 rolls | 12 cm | Blue | 4038 |
| | | Black | 4038Z |



ORFICAST® MORE 15 CM

ORFICAST® MORE 15 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|--------|---------|
| 1 roll | 15 cm | Blue | 4039 |
| | | Black | 4039Z |
| | | Orange | 4039OR |
| 6 rolls | 15 cm | Blue | 4040 |
| | | Black | 4040Z |
| | | Orange | 4040OR |



ORFICAST® MORE 30 CM

ORFICAST® MORE 30 CM - LENGTH 3 M

| Quantity | Width | Color | Art. N° |
|----------|-------|-------|---------|
| 1 roll | 30 cm | Blue | 4042 |




ORFIT® CLASSIC

ORFIT® CLASSIC SOFT - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 1.6 mm (1/16") | non perfo | 8332.S01 | |
| | micro perfo | 8332.S02 | |
| 2.0 mm (1/12") | non perfo | 8333.S01 | |
| | micro perfo | 8333.S02 | |
| | micro plus perfo | 8333.S02+ | |
| | maxi perfo | 8333.S03 | |
| | mini perfo | 8333.S04 | |
| 2.5 mm (3/32") | micro perfo | 8338.S02 | |
| 3.2 mm (1/8") | non perfo | 8334.S01 | 8354.S01 |
| | maxi perfo | 8334.S03 | 8354.S03 |
| | mini perfo | 8334.S04 | 8354.S04 |
| 4.2 mm (3/16") | non perfo | | 8355.S01 |
| | mini perfo | | 8355.S04 |

ORFIT® CLASSIC STIFF - SHEET MATERIAL


| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 3.2 mm (1/8") | non perfo | 8334.ST1 | 8354.ST1 |
| | mini perfo | 8334.ST4 | 8354.ST4 |
| 4.2 mm (3/16") | non perfo | | 8355.ST1 |
| | mini perfo | | 8355.ST4 |

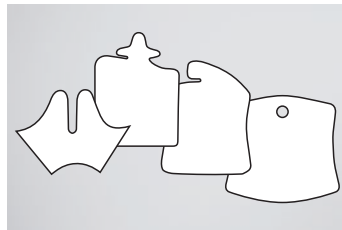

ORFIT® CLASSIC PRECUTS

ORFIT® CLASSIC PRECUTS (packed per 2 pieces)

| | | Size | Art. N° |
|----------------------------------|----------------------------------|---------|---------|
| GAUNTLET THUMB POST ORTHOSIS | | | |
| Orfit Classic | 1.6 mm {1/16"} micro perfo, soft | small | 35810 |
| | | medium | 35811 |
| | | large | 35812 |
| COMPLETE THUMB ORTHOSIS | | | |
| Orfit Classic | 2.0 mm {1/12"} micro perfo, soft | small | 35900 |
| | | medium | 35901 |
| | | large | 35902 |
| WRIST + THUMB ORTHOSIS | | | |
| Orfit Classic | 2.0 mm {1/12"} micro perfo, soft | small | 35814 |
| | | medium | 35815 |
| | | large | 35816 |
| GAUNTLET IMMOBILIZATION ORTHOSIS | | | |
| Orfit Classic | 2.0 mm {1/12"} micro perfo, soft | small | 35830 |
| | | medium | 35831 |
| | | large | 35832 |
| DORSAL WRIST COCK-UP ORTHOSIS | | | |
| Orfit Classic | 3.2 mm {1/8"} non perfo, stiff | x-small | 35818 |
| | | small | 35820 |
| | | medium | 35821 |
| | | large | 35822 |

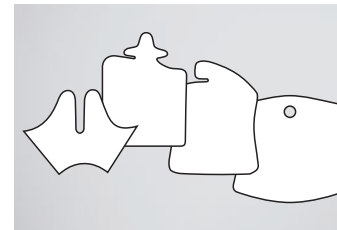

ORFIT® CLASSIC PRECUTS

ORFIT® CLASSIC PRECUTS (packed per 2 pieces)

| | | Size | Art. N° |
|---------------------------------------|---------------------------------|--------|---------|
| INTRINSIC RESTING HAND ORTHOSIS | | | |
| Orfit Classic | 3.2 mm (1/8") non perfo, stiff | small | 35850 |
| | | medium | 35851 |
| | | large | 35852 |
| ANTI-SPASTIC ORTHOSIS + THUMB PIECE | | | |
| Orfit Classic | 3.2 mm (1/8") non perfo, stiff | small | 35870 |
| | | medium | 35871 |
| | | large | 35872 |
| RHEUMATOID ARTHRITIS RESTING ORTHOSIS | | | |
| Orfit Classic | 3.2 mm (1/8") mini perfo, stiff | small | 35840 |
| | | medium | 35841 |
| | | large | 35842 |
| ULNAR DRIFT ORTHOSIS | | | |
| Orfit Classic | 3.2 mm (1/8") non perfo, stiff | left | 35817KL |
| | | right | 35817KR |



ORFIT® NS

ORFIT® NS SOFT - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 1.6 mm (1/16") | non perfo | 8332.S01/NS | |
| | micro perfo | 8332.S02/NS | |
| 2.0 mm (1/12") | non perfo | 8333.S01/NS | |
| | micro perfo | 8333.S02/NS | |
| | micro plus perfo | 8333.S02+/NS | |
| | maxi perfo | 8333.S03/NS | |
| | mini perfo | 8333.S04/NS | |
| 2.5 mm (3/32") | micro perfo | 8338.S02/NS | |
| 3.2 mm (1/8") | non perfo | 8334.S01/NS | 8354.S01/NS |
| | maxi perfo | 8334.S03/NS | 8354.S03/NS |
| | mini perfo | 8334.S04/NS | 8354.S04/NS |
| 4.2 mm (3/16") | non perfo | | 8355.S01/NS |
| | mini perfo | | 8355.S04/NS |

ORFIT® NS STIFF - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 3.2 mm (1/8") | non perfo | 8334.ST1/NS | 8354.ST1/NS |
| | mini perfo | 8334.ST4/NS | 8354.ST4/NS |
| 4.2 mm (3/16") | non perfo | | 8355.ST1/NS |
| | mini perfo | | 8355.ST4/NS |



ORFIT® COLORS NS

ORFIT® COLORS NS - SHEET MATERIAL


| Thickness (mm-inch) | Perforation types | Color | 45 cm x 60 cm (17 11/16" x 23 5/8") |
|------------------------|----------------------|----------------|--|
| 2.0 mm (1/12") | non perfo | Violet | 8133PU.1/NS |
| | | Hot Green | 8133GR.1/NS |
| | | Dominant Black | 8133Z.1/NS |
| | | Dynamic Red | 8133R.1/NS |
| | | Ocean Blue | 8133B.1/NS |
| | | Bright Pink | 8133FP.1/NS |
| | | Sonic Silver | 8133MG.1/NS |
| | | Atomic Blue | 8133MB.1/NS |
| | | Gold | 8133GO.1/NS |
| 2.0 mm (1/12") | micro perfo | Violet | 8133PU.2/NS |
| | | Hot Green | 8133GR.2/NS |
| | | Dominant Black | 8133Z.2/NS |
| | | Dynamic Red | 8133R.2/NS |
| | | Ocean Blue | 8133B.2/NS |
| | | Bright Pink | 8133FP.2/NS |
| | | Sonic Silver | 8133MG.2/NS |
| | | Atomic Blue | 8133MB.2/NS |
| | | Gold | 8133GO.2/NS |



ORFIT® COLORS NS

ORFIT® COLORS NS - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | Color | 45 cm x 60 cm (17 11/16" x 23 5/8") |
|------------------------|----------------------|--------------|--|
| 3.4 mm (1/8") | non perfo | Sonic Silver | 8134MG.1/NS |
| | | Atomic Blue | 8134MB.1/NS |
| | | Gold | 8134GO.1/NS |
| 3.4 mm (1/8") | mini perfo | Sonic Silver | 8134MG.4/NS |
| | | Atomic Blue | 8134MB.4/NS |
| | | Gold | 8134GO.4/NS |

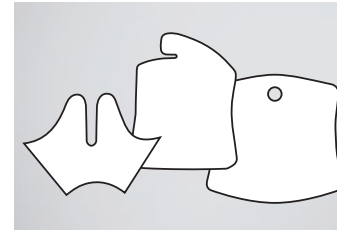

ORFIT® COLORS NS PRECUTS

ORFIT® COLORS NS PRECUTS (packed per 2 pieces)

| Color | Small | Medium | Large |
|---|------------|------------|------------|
| GAUNTLET THUMB POST ORTHOSIS | | | |
| Ocean Blue | 35810B/NS | 35811B/NS | 35812B/NS |
| Atomic Blue | 35810MB/NS | 35811MB/NS | 35812MB/NS |
| Sonic Silver | 35810MG/NS | 35811MG/NS | 35812MG/NS |
| Dominant Black | 35810Z/NS | 35811Z/NS | 35812Z/NS |
| WRIST + THUMB ORTHOSIS | | | |
| Ocean Blue | 35814B/NS | 35815B/NS | 35816B/NS |
| Atomic Blue | 35814MB/NS | 35815MB/NS | 35816MB/NS |
| Sonic Silver | 35814MG/NS | 35815MG/NS | 35816MG/NS |
| Dominant Black | 35814Z/NS | 35815Z/NS | 35816Z/NS |
| GAUNTLET IMMOBILIZATION ORTHOSIS | | | |
| Ocean Blue | 35830B/NS | 35831B/NS | 35832B/NS |
| Atomic Blue | 35830MB/NS | 35831MB/NS | 35832MB/NS |
| Sonic Silver | 35830MG/NS | 35831MG/NS | 35832MG/NS |
| Dominant Black | 35830Z/NS | 35831Z/NS | 35832Z/NS |



ORFIT® STRIPS

ORFIT® STRIPS - SHEET MATERIAL - THIN

| Thickness | Pieces / Color | 45 cm x 0.5 cm | 45 cm x 2 cm |
|----------------|-----------------------------|----------------|--------------|
| 2.0 mm (1/12") | 10 x Sonic Silver | 8113MG.1/NS | |
| | 10 x Gold | 8113GO.1/NS | |
| | 5 x Sonic Silver - 5 x Gold | 8113GOMG.1/NS | |

ORFIT® STRIPS - SHEET MATERIAL - WIDE

| Thickness | Pieces / Color | 45 cm x 0.5 cm | 45 cm x 2 cm |
|---------------|----------------------------|----------------|---------------|
| 3.4 mm (1/8") | 10 x Gold | | 8124GO.1/NS |
| | 10 x Atomic Blue | | 8124MB.1/NS |
| | 5 x Gold - 5 x Atomic Blue | | 8124MBGO.1/NS |



ORFILIGHT® NS ORFILIGHT® BLACK NS ORFILIGHT® ATOMIC BLUE NS

ORFILIGHT® NS - SHEET MATERIAL


| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") |
|---------------------|-------------------|-------------------------------------|
| 1.6 mm (1/16") | micro perfo | 8332.2/L |
| 2.5 mm (3/32") | micro perfo | 8338.2/L |
| 3.2 mm (1/8") | non perfo | 8334.1/L |
| | mini perfo | 8334.4/L |

ORFILIGHT® BLACK NS - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | Color | 45 cm x 60 cm (17 11/16" x 23 5/8") |
|---------------------|-------------------|-------|-------------------------------------|
| 1.6 mm (1/16") | micro perfo | Black | 8332Z.2/L/NS |
| 2.5 mm (3/32") | micro perfo | Black | 8338Z.2/L/NS |
| 3.2 mm (1/8") | non perfo | Black | 8334Z.1/L/NS |
| | mini perfo | Black | 8334Z.4/L/NS |

ORFILIGHT® ATOMIC BLUE NS - SHEET MATERIAL


| Thickness (mm-inch) | Perforation types | Color | 45 cm x 60 cm (17 11/16" x 23 5/8") |
|---------------------|-------------------|-------------|-------------------------------------|
| 1.6 mm (1/16") | micro perfo | Atomic Blue | 8332MB.2/L/NS |
| 2.5 mm (3/32") | micro perfo | Atomic Blue | 8338MB.2/L/NS |
| 3.2 mm (1/8") | non perfo | Atomic Blue | 8334MB.1/L/NS |
| | mini perfo | Atomic Blue | 8334MB.4/L/NS |

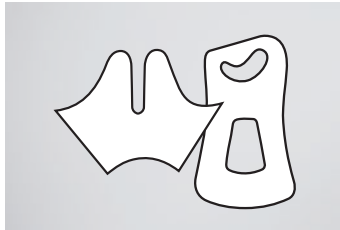

ORFILIGHT® NS PRECUTS

ORFILIGHT® NS PRECUTS (packed per 2 pieces)

| | | Size | Art. N° |
|-------------------------------|----------------------------|--------|---------|
| GAUNTLET THUMB POST ORTHOSIS | | | |
| Orfilight NS | 1.6 mm [1/16"] micro perfo | small | 35810/L |
| | | medium | 35811/L |
| | | large | 35812/L |
| DORSAL WRIST COCK-UP ORTHOSIS | | | |
| Orfilight NS | 3.2 mm [1/8"] non perfo | small | 35820/L |
| | | medium | 35821/L |
| | | large | 35822/L |

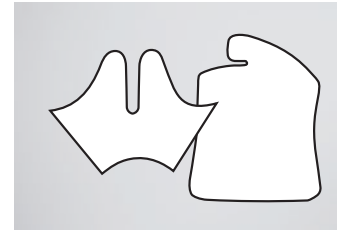

ORFILIGHT® BLACK NS PRECUTS

ORFILIGHT® BLACK NS PRECUTS (packed per 2 pieces)

| | | Size | Art. N° |
|------------------------------|----------------------------|--------|-------------|
| GAUNTLET THUMB POST ORTHOSIS | | | |
| Orfilight Black NS | 1.6 mm [1/16"] micro perfo | small | 35810Z/L/NS |
| | | medium | 35811Z/L/NS |
| | | large | 35812Z/L/NS |
| WRIST + THUMB ORTHOSIS | | | |
| Orfilight Black NS | 2.5 mm [3/32"] micro perfo | small | 35814Z/L/NS |
| | | medium | 35815Z/L/NS |
| | | large | 35816Z/L/NS |



ORFIGUARD

ORFIGUARD

| Thickness | Color | one-size fits all |
|----------------|----------------|-------------------|
| 2.0 mm (1/12") | Dominant Black | 35903Z/NS |



AQUAFIT™ NS (ORFIT® NATURAL NS)

AQUAFIT™ NS (ORFIT® NATURAL NS) - SOFT - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 1.6 mm (1/16") | non perfo | 8632N.1/NS | |
| | micro perfo | 8632N.2/NS | |
| 2.4 mm (3/32") | non perfo | 8639N.1/NS | |
| | micro plus perfo | 8639N.2+/NS | |
| 3.2 mm (1/8") | non perfo | 8634N.1/NS | |
| | macro perfo | 8634N.6/NS | |
| | maxi perfo | 8634N.3/NS | |
| | non perfo | | 8654N.1/NS |
| | macro perfo | | 8654N.6/NS |

AQUAFIT™ NS (ORFIT® NATURAL NS) - STIFF - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 3.2 mm (1/8") | non perfo | 8634N.ST1/NS | 8654N.ST1/NS |
| | macro perfo | 8634N.ST6/NS | 8654N.ST6/NS |



TECNOFIT™

TECNOFIT™ - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 2.5 mm (3/32") | non perfo | 8738.1 | |
| | maxi perfo | 8738.3 | |
| | mini perfo | 8738.4 | |
| 3.2 mm (1/8") | non perfo | 8734.1 | 8754.1 |
| | mini perfo | 8734.4 | 8754.4 |



ORFIZIP® LIGHT NS

ORFIZIP® LIGHT NS IVORY

| Thickness (mm-inch) | Perforation types | Size | Art. N° |
|------------------------|----------------------|--------|------------------|
| 2.5 mm (3/32") | micro perfo | small | 39108N/25MI/L/NS |
| | | medium | 39109N/25MI/L/NS |
| | | large | 39110N/25MI/L/NS |
| 3.2 mm (1/8") | | small | 39111N/32MI/L/NS |
| | | medium | 39112N/32MI/L/NS |
| | | large | 39113N/32MI/L/NS |

ORFIZIP® LIGHT NS ATOMIC BLUE

| Thickness (mm-inch) | Perforation types | Size | Art. N° |
|------------------------|----------------------|--------|-------------------|
| 2.5 mm (3/32") | micro perfo | small | 39108MB/25MI/L/NS |
| | | medium | 39109MB/25MI/L/NS |
| | | large | 39110MB/25MI/L/NS |
| 3.2 mm (1/8") | micro perfo | small | 39111MB/32MI/L/NS |
| | | medium | 39112MB/32MI/L/NS |
| | | large | 39113MB/32MI/L/NS |



MULTIFIT® NS

MULTIFIT® NS PEARL WHITE - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 3.2 mm (1/8") | non perfo | 7134.1/NS | |
| | non perfo | | 7154.1/NS |
| | macro perfo | 7134.6/NS | |
| | macro perfo | | 7154.6/NS |

MULTIFIT® NS BLACK - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 3.2 mm (1/8") | non perfo | 7134Z.1/NS | |
| | non perfo | | 7154Z.1/NS |
| | macro perfo | 7134Z.6/NS | |
| | macro perfo | | 7154Z.6/NS |




ORFIT® ECO ORFIT® ECO BLACK NS

ORFIT® ECO - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 2.4 mm (3/32") | non perfo | 8938N.1 | |
| | mini perfo | 8938N.4 | |
| 3.2 mm (1/8") | non perfo | 8934N.1 | 8954N.1 |
| | mini perfo | 8934N.4 | 8954N.4 |

ORFIT® ECO BLACK NS - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 45 cm x 60 cm (17 11/16" x 23 5/8") | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--|--------------------------------------|
| 2.4 mm (3/32") | non perfo | 8938Z.1/NS | |
| | mini perfo | 8938Z.4/NS | |
| 3.2 mm (1/8") | non perfo | 8934Z.1/NS | 8954Z.1/NS |
| | mini perfo | 8934Z.4/NS | 8954Z.4/NS |

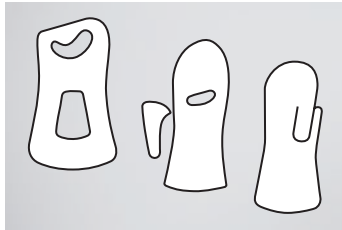

ORFIT® ECO PRECUTS

ORFIT® ECO PRECUTS (packed per 2 pieces)

| | | Size | Art. N° |
|--------------------------------------|--------------------------|--------|---------------|
| DORSAL WRIST COCK-UP ORTHOSIS | | | |
| Orfit Eco | 3.2 mm (1/8") non perfo | small | 35820/ECO |
| | | medium | 35821/ECO |
| | | large | 35822/ECO |
| ANTI-SPASTIC ORTHOSIS + THUMB PIECE | | | |
| Orfit Eco Black NS | 3.2 mm (1/8") mini perfo | small | 35870Z/ECO/NS |
| | | medium | 35871Z/ECO/NS |
| | | large | 35872Z/ECO/NS |
| INTRINSIC ANTI-SPASTIC HAND ORTHOSIS | | | |
| Orfit Eco Black NS | 3.2 mm (1/8") mini perfo | small | 35910Z/ECO/NS |
| | | medium | 35911Z/ECO/NS |
| | | large | 35912Z/ECO/NS |



ORFIBRACE™ NS

ORFIBRACE™ NS - SHEET MATERIAL

| Thickness (mm-inch) | Perforation types | 60 cm x 90 cm (23 5/8" x 35 3/8") |
|------------------------|----------------------|--------------------------------------|
| 3.0 mm (1/8") | non perfo | 8854.1/NS |
| | macro perfo | 8854.6/NS |



ALLFIT™

ALLFIT™

| Quantity | Art. N° |
|--------------------------|---------|
| Box with 6 bags of 500 g | 703011 |



DYNASYST™

DYNASYST™ ACCESSORIES

| | Art. N° |
|---|---------|
| ORFITUBE | |
| Orfitube™ kit | 35311 |
| Orfitube™ adjuster kit | 35305 |
| Orfitube™ ends | 35318 |
| Orfitubes™ + Bending wires | 35310 |
| Single Orfitube™ adjuster - 5 pcs. + Allen wrench | 35304 |
| Orfitube™ bending tool | 35308 |
| COILS AND SPRINGS | |
| Knuckle bender coil spring | 35801 |
| Finger coil spring | 35800 |
| Finger extension assist | 35802 |
| Wrist extension coil spring | 35804 |
| Radial wrist extension splint (small) | 35865K |
| Radial wrist extension splint (medium) | 35866K |
| Spring wire | 35805 |



DYNASYST™

DYNASYST™ ACCESSORIES

| | Art. N° |
|---|---------|
| ISOFORCE | |
| Isoforce (small) | 35322 |
| Isoforce (medium) | 35323 |
| Isoforce (large) | 35324 |
| THREADS | |
| Polyester monofilament thread | 35315 |
| Elastic thread (max. resistance 200 g.) | 35312 |
| Elastic thread (max. resistance 400 g.) | 35313 |
| Elastic thread (max. resistance 600 g.) | 35314 |

WORKING TOOLS

| | Art. N° |
|-------------------------------------|----------|
| Sharp Scissors - Curved (14 cm) | 35061 |
| Strong scissors | 35065 |
| Gingher shears | 35066 |
| Bandage scissors | 35067 |
| Spring coil pliers | 35140 |
| Chrome plated eyelet setter | 35134 |
| 200 eyelets of 4 mm (1/6" diameter) | 35137 |
| Deburring knife | 35307 |
| 10 spare blades with straight cut | 35307B10 |
| Spare blades with oblique cut | 35307B20 |
| Dynamometer | 35105 |
| Yellow marker | 35003G |
| Refill yellow marker | 35003G/1 |
| Black marker | 35003Z |
| Refill black marker | 35003Z/1 |
| White marker | 35003W |
| Refill white marker | 35003W/1 |

PADDING AND LINING

| | Art. N° |
|----------------------------------|-------------|
| Luxofoam non perforated | 35110 |
| Luxofoam perforated | 35112 |
| Orfilastic | 35111 |
| Orfilastic Black | 35117 |
| Edging strip (natural color) | 360201N.ST1 |
| Soft velvet edging strip (black) | 35116 |

HOOK AND LOOP STRAPS - 25 MM

| STANDARD LOOP STRAP | Length | Color | Art. N° |
|--------------------------------------|----------------|-------------|---------|
| Standard loop strap - 25 mm (1") | 25 m (98 2/5") | white | 35400L |
| | | beige | 35411L |
| | | blue | 35401L |
| | | black | 35404L |
| | | fluo yellow | 35402L |
| | | fluo pink | 35403L |
| SELF-ADHESIVE LOOP TAPE | Length | Color | Art. N° |
| Self-adhesive loop tape - 25 mm (1") | 25 m (98 2/5") | white | 35405L |
| | | beige | 35412L |
| STANDARD HOOK STRAP | Length | Color | Art. N° |
| Standard hook strap - 25 mm (1") | 25 m (98 2/5") | white | 35400H |
| | | beige | 35411H |
| | | black | 35417H |
| SELF-ADHESIVE HOOK TAPE | Length | Color | Art. N° |
| Self-adhesive hook tape - 25 mm (1") | 25 m (98 2/5") | white | 35405H |
| | | beige | 35412H |
| | | blue | 35406H |
| | | black | 35409H |
| | | fluo yellow | 35407H |
| | | fluo pink | 35408H |

HOOK AND LOOP STRAPS - 50 MM

| STANDARD LOOP STRAP | Length | Color | Art. N° |
|--------------------------------------|----------------|-------|---------|
| Standard loop strap - 50 mm (2") | 25 m (98 2/5") | white | 35413L |
| | | beige | 35415L |
| | | black | 35441L |
| SELF-ADHESIVE LOOP TAPE | | Color | Art. N° |
| Self-Adhesive loop tape - 50 mm (2") | 25 m (98 2/5") | white | 35414L |
| | | beige | 35416L |
| STANDARD HOOK STRAP | | Color | Art. N° |
| Standard hook strap - 50 mm (2") | 25 m (98 2/5") | white | 35413H |
| | | beige | 35415H |
| SELF-ADHESIVE HOOK TAPE | | Color | Art. N° |
| Self-adhesive hook tape - 50 mm (2") | 25 m (98 2/5") | white | 35414H |
| | | beige | 35416H |
| | | black | 35442H |

DOUBLE-SIDED HOOK OR LOOP STRAP

| | Width | Length | Color | Art. N° |
|--------------------------------|------------|------------|-------|---------|
| Double-sided velvet loop strap | 25 mm (1") | 25 m (82') | black | 35444L |
| Double-sided micro hook strap | 25 mm (1") | 25 m (82') | black | 35444H |

ELASTIC LOOP TAPE

| Width | Length | Color | Art. N° |
|------------|----------------|-------|---------|
| 25 mm (1") | 25 m (98 2/5") | white | 35437L |
| | | beige | 35439L |
| | | black | 35432L |
| 50 mm (2") | 25 m (98 2/5") | white | 35438L |
| | | beige | 35440L |
| | | black | 35443L |

SPLINT STRAP

| Width | Length | Color | Art. N° |
|------------|------------------|-------|---------|
| 25 mm (1") | 250 mm (10") | beige | 35450 |
| 25 mm (1") | 400 mm (15 3/4") | beige | 35455 |
| 25 mm (1") | 500 mm (20") | beige | 35460 |

HEATING DEVICES

| | Art. N° |
|-----------------------------|-----------------------|
| Water Bath | 32502/230EU |
| Large Water Bath (Suspan 1) | 35097/115/1400 (115V) |
| Large Water Bath (Suspan 1) | 35097/230/1400 (220V) |
| Strava Splint Bath | 35123/230EU |
| Dry Heater | 32503/230EU |
| Dry Heater | 32503/230UK |
| Dry Heater | 32503/230CH |
| Heat Gun | 35353 |

| | Properties | | | | Thickness | | | | | | Available in precuts | Transparent when activated |
|---------------------------------------|------------|-------------|----------|-----------|-----------|--------|--------|--------|--------|--------|-------------------------|-------------------------------|
| | Elastic | Low stretch | Flexible | On a roll | 1.6 mm | 2.0 mm | 2.4 mm | 2.5 mm | 3.2 mm | 4.2 mm | | |
| Orficast® | | | | ■ | | | | | | | | |
| Orfit® Classic Soft | ■ | | | | ■ | ■ | | ■ | ■ | ■ | ■ | ■ |
| Orfit® Classic Stiff | ■ | | | | | | | | ■ | ■ | | ■ |
| Orfit® NS Soft | ■ | | | | ■ | ■ | | ■ | ■ | ■ | | ■ |
| Orfit® NS Stiff | ■ | | | | | | | | ■ | ■ | | ■ |
| Orfit® Colors NS | ■ | | | | | ■ | | | ■ | | ■ | |
| Orfilight® NS | ■ | | | | ■ | | | ■ | ■ | | ■ | |
| Orfilight® Black NS | ■ | | | | ■ | | | ■ | ■ | | ■ | |
| Orfilight® Atomic Blue NS | ■ | | | | ■ | | | ■ | ■ | | | |
| Aquafit™ NS Soft (Orfit® Natural NS) | | ■ | | | | ■ | | | ■ | ■ | | ■ |
| Aquafit™ NS Stiff (Orfit® Natural NS) | | ■ | | | | | | | | ■ | | ■ |
| Tecnofit™ | ■ | | | | | | | ■ | ■ | | | ■ |
| Orfizip® Light NS | ■ | | | | | | | ■ | ■ | | ■ | |
| Multifit® NS | | ■ | | | | | ■ | | ■ | | | |
| Orfit® Eco | | ■ | | | | | | ■ | ■ | | ■ | |
| Orfit® Eco Black NS | | ■ | | | | | | ■ | ■ | | ■ | |
| Orfibrace™ NS | | ■ | | | | | | | ■ | | ■ | |

| Resistance to stretch | | | Memory | | | | Drape | | | | | | |
|-----------------------|--------|------|--------|----------|------|------|-------|----------|------|-------------------|------------------------|-----------------|-----------------|
| Low | Medium | High | Low | Moderate | High | Full | Low | Moderate | High | Non-stick Coating | Activation Temperature | Working time | Hardening time |
| | ■ | | | ■ | | | | | ■ | | 65°C-149°F | 1.0-1.5 minutes | 3-3.5 minutes |
| | ■ | | | | | ■ | | ■ | | | 65°C-149°F | 1- 2 minutes | 1-6 minutes |
| | | ■ | | | | ■ | ■ | | | | 65°C-149°F | 2-3 minutes | 6-10 minutes |
| | ■ | | | | | ■ | | ■ | | ■ | 65°C-149°F | 1- 2 minutes | 2-6 minutes |
| | | ■ | | | | ■ | ■ | | | ■ | 65°C-149°F | 2-3 minutes | 6-9 minutes |
| | ■ | | | | | ■ | | ■ | | ■ | 65°C-149°F | 1-3.5 minutes | 3.5-9 minutes |
| | ■ | | | | ■ | ■ | | ■ | | ■ | 65°C-149°F | 1-2.5 minutes | 2-4 minutes |
| | ■ | | | | ■ | | | ■ | | ■ | 65°C-149°F | 1-2.5 minutes | 2-4 minutes |
| | ■ | | | | ■ | | | ■ | | ■ | 65°C-149°F | 1-2.5 minutes | 2-4 minutes |
| | ■ | | | | | ■ | | ■ | | ■ | 65°C-149°F | 1-2.5minutes | 3-6 minutes |
| | | ■ | | | | ■ | ■ | | | ■ | 65°C-149°F | 1.5-2.5 minutes | 8-9.5 minutes |
| | | ■ | | | | ■ | ■ | | | | 65°C-149°F | 1-3 minutes | 2.5-5.5 minutes |
| | ■ | | | ■ | | | | ■ | | ■ | 65°C-149°F | 1-1.5 minutes | 2-4 minutes |
| | ■ | | ■ | | | | | ■ | | ■ | 65°C-149°F | 2 minutes | 2 minutes |
| | | ■ | | ■ | | | | ■ | | | 65°C-149°F | 1-3 minutes | 4-6 minutes |
| | | ■ | | ■ | | | | ■ | | ■ | 65°C-149°F | 1-3 minutes | 4-6 minutes |
| | | ■ | | ■ | | | | ■ | | ■ | 65°C-149°F | 1-2 minutes | 5-6 minutes |

Our Commitment to Education and Inspiration

The Orfit Academy

The Orfit Academy is a virtual school in which we offer multiple orthotic fabrication courses on various themes, covering all aspects of hand rehabilitation. The courses offer a detailed look at the anatomy, the clinical diagnoses and conditions, multiple fabrication videos and orthotic designs, case studies, and evidence-based resources.

Go to orfit.com/physical-rehabilitation/academy and enroll today for free access to multiple courses that will inspire you to be the best clinician possible for your patients

Daily inspiration

Follow the social media channel of your choice for daily inspiration, splinting tips and tricks.



[@orfitindustries](https://www.instagram.com/orfitindustries)



[Orfit Splinting & Physical Rehabilitation](https://www.linkedin.com/company/orfit-splinting-physical-rehabilitation)

Connect with other therapists

The community of Orfit users spans across the world and they love to share their most interesting patient cases and creative orthotic fabrication ideas in the Orfit Splinting Group on Facebook. Join us now, learn from each other and support your colleagues.



facebook.com/groups/Orfit.splinting

In-person or online workshops



Contact Debby Schwartz at debby.schwartz@orfit.com to see how to set up a training at your facility and earn valuable continuing education credit. Select from a variety of immobilization and/or mobilization orthoses to create a custom workshop that meets your clinical requirements.




Orfit works closely together with its distributors to offer orthotic fabrication workshops that enhance your skills and introduce new innovative ideas in the world of orthotic fabrication. Contact us to inform about workshops in your area. Go to orfit.com/contact to submit a workshop request.



Art. N° 21021E - THE SPLINTING GUIDE

Compiled by Paul Van Lede OT, MS.

- > Static hand orthoses.
- > Dynamic hand orthoses.

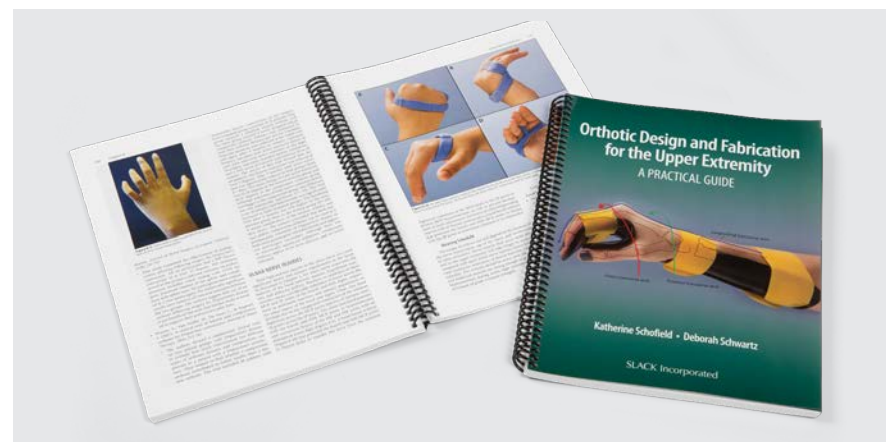

Art. N° 52013 - ORTHOTIC DESIGN AND FABRICATION FOR THE UPPER EXTREMITY

A Practical Guide by Orfit Product Specialist Debby Schwartz and Drs. Katherine Schofield covers foundational knowledge for the use of orthoses as an intervention in occupational therapy practice. The combination of illustrations, text and instructional video material makes Orthotic Design and Fabrication for the Upper Extremity a valuable resource for occupational therapy students, new graduates, and novice clinicians.

Our Commitment to Innovation and Quality

Quality Assurance Management

Orfit Industries permanently invests in research and development to provide innovative and sustainable products with the highest quality level. We work in close collaboration with leading rehabilitation centers throughout the world, developing solutions that improve your orthotic fabrication activities and the clinical outcomes for your patients.

Orthotic applications in physical rehabilitation require accuracy and durability. The same goes for our production process.

Our commitment to providing excellent quality products is backed up by compliance with ISO 13485 standards, as well as with European CE, MDR and other regulations.



We permanently and systematically update our quality control as proven by our ISO 13485:2016 certification. The certification demonstrates our ability to provide medical devices that consistently meet customer and regulatory requirements.




The European Union Medical Device Regulation entered into application on 26 May 2021. It is a set of regulations that governs the production and distribution of medical devices in Europe. A suite of MDR-compliant regulatory systems, processes and documents continually monitor the safety and performance of our products.



FDA stands for the U.S. Food and Drug Administration, an agency of the United States Department of Health and Human services. The FDA legislation protects public health through the regulation and supervision of food and medical products. Our products are in compliance with the Federal Food, Drug and Cosmetic Act Code of Federal Regulations Title 21 part 808/812/820.



Products that are brought onto the European market must meet the essential requirements of the European product CE label certification regulation. With this label, we declare that all products made by Orfit comply with the essential requirements of the European health, safety and environmental protection legislation.

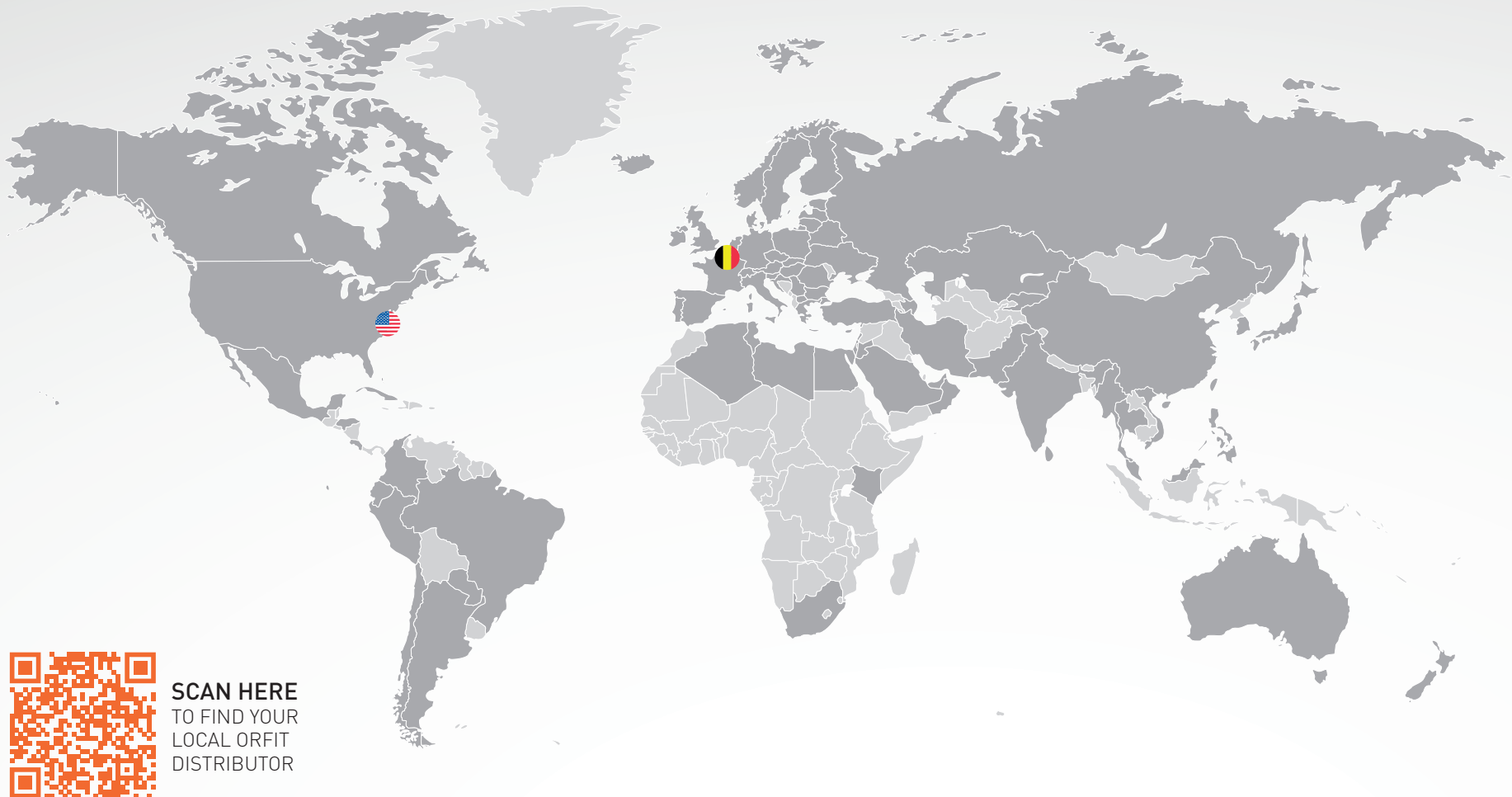

SCAN HERE
TO FIND YOUR
LOCAL ORFIT
DISTRIBUTOR

orfit.com/contact

Offices in 2 countries

Sales in +90 countries

Network of 89 distributors



www.orfit.com



ORFIT INDUSTRIES

Vosveld 9A | B-2110 Wijnegem | Belgium
T (+32) (0)3 326 20 26
welcome@orfit.com

ORFIT INDUSTRIES AMERICA

810 Ford Drive | Norfolk, VA 23523
T 516-935-8500
welcome@orfit.com

10-11-2022
Art. N° 520001