


LMB Spring Finger Extension Assist

Indicated for PIP joint limitation of 45° or less. Approximately 8oz. (227g) of force is applied when the PIP joint is in 15° of flexion. The contoured center pad distributes pressure evenly on the PIP joint. Bilateral wires may be bent to adjust force and to accommodate changes in edema. To size, measure from the distal palmar crease to the distal finger crease (finger extended).

+ Treat a variety of diagnoses involving PIP joint tightness

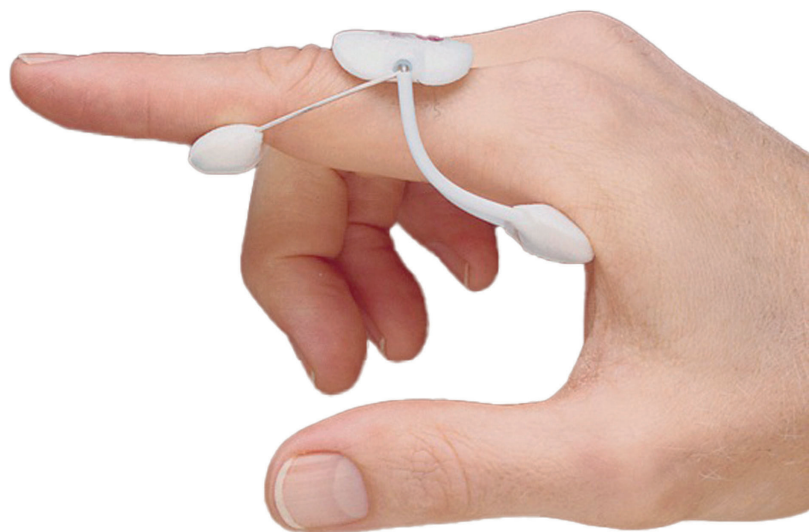

Size	Code	DPC to DIP Length
XS	501A	5.7cm
S	501AA	6.7cm
M	501B	7.6cm
L	501C	8.9cm
XL	501D	10cm